

HPLC Application

ID No.: 25472

Aloe Emodin, Eugenol Methyl Ether, and Aloins (A and B) from Aloe Standards Screen on Kinetex 2.6u PS C18 50x

Column: Kinetex® 2.6 µm PS C18 100 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4780-E0

Elution Type: Gradient

Eluent A: 0.1% Formic Acid in Water

Eluent B: 0.1% Formic Acid in MeOH

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25
	2	0.5	75	25
	3	5	20	80
	4	7	20	80
	5	7.1	75	25
	6	9	75	25


Flow Rate: 1.2 mL/min

Col. Temp.: 30 °C

Detection: LC/UV (DAD, PDA) @ 220.000000000 nm (nanometers) (30 °C)
25472



Products used in this application:

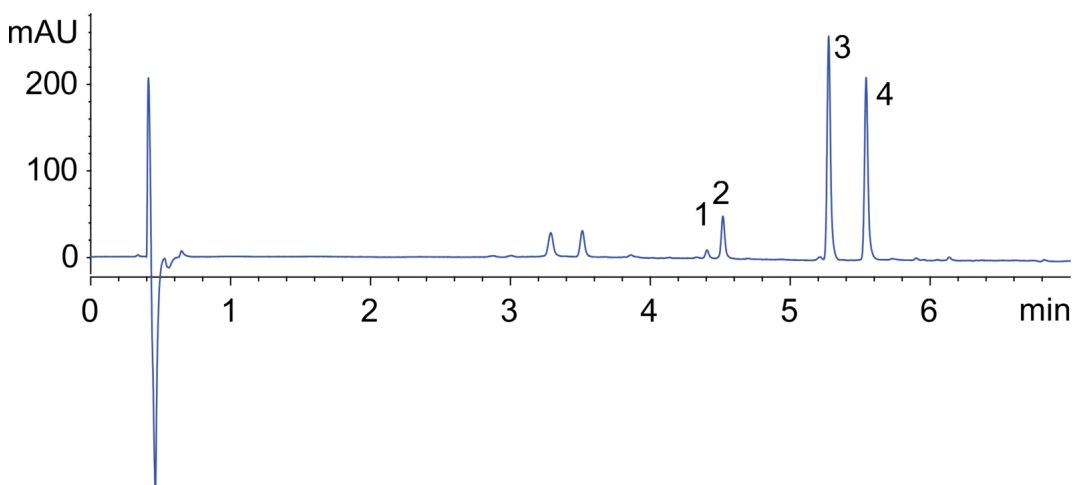

ANALYTES:

- 1 Aloin A
Retention Time: 4.405 min
- 2 Aloin B
Retention Time: 4.519 min
- 3 Methyl Eugenol
Retention Time: 5.275 min
- 4 Aloe Emodin
Retention Time: 5.543 min

