

Analysis of Uracil and metabolites on Kinetex 2.6 µm PS C18 150x3.0 Column

Column: Kinetex® 2.6 µm PS C18 100 Å, LC Column 150 x 3 mm, Ea

Dimensions: 150 x 3 mm ID

Order No: 00F-4780-Y0

Elution Type: Gradient

Eluent A: 0.1% Formic acid in Water

Eluent B: Methanol


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	93	7
	2	6	93	7
	3	6.01	50	50
	4	7	50	50
	5	7.01	93	7
	6	12	93	7

Flow Rate: 0.2 mL/min


Col. Temp.: 15 °C

Detection: LC/MS/MS @ 0.0000000000 m/z (mass per charge) (15 °C)

Analyst Note: Analyte diluted with initial mobile phase

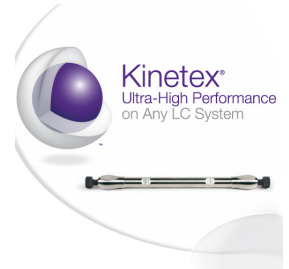

Products used in this application:

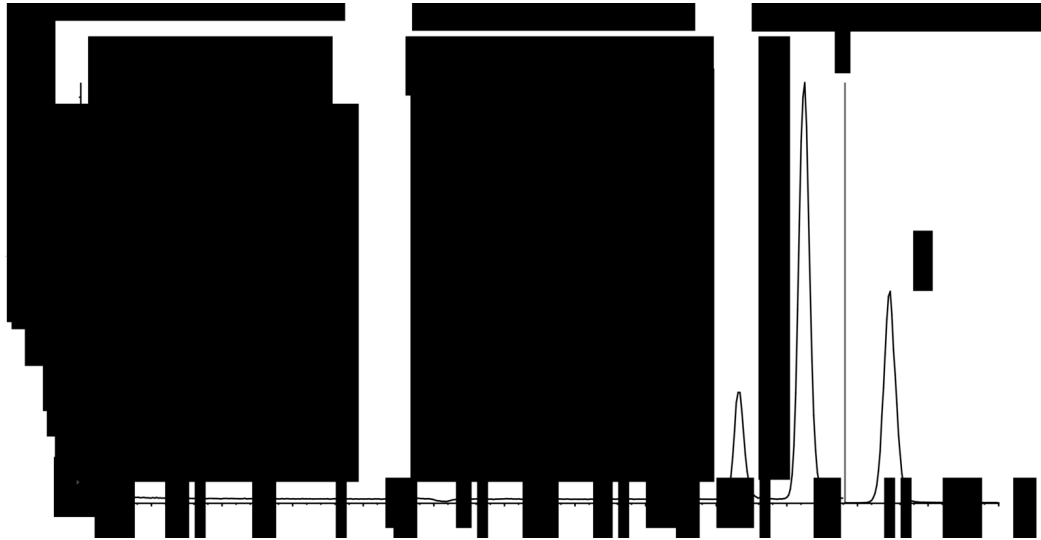

ANALYTES:

1 5,6-dihydrouracil (UH2)

Retention Time: 4.6 min

2 Uracil

Retention Time: 5.1 min

3 5-fluorouracil

Retention Time: 5.7 min

