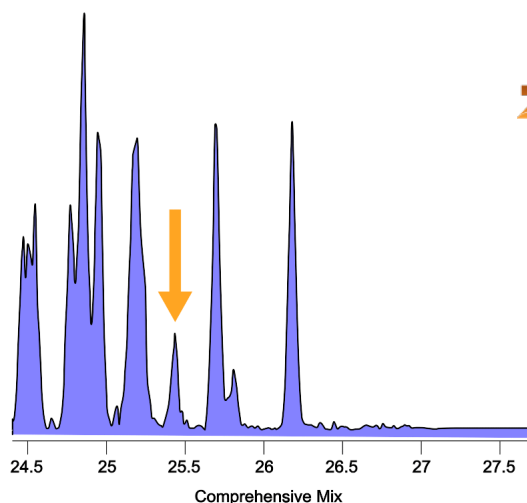
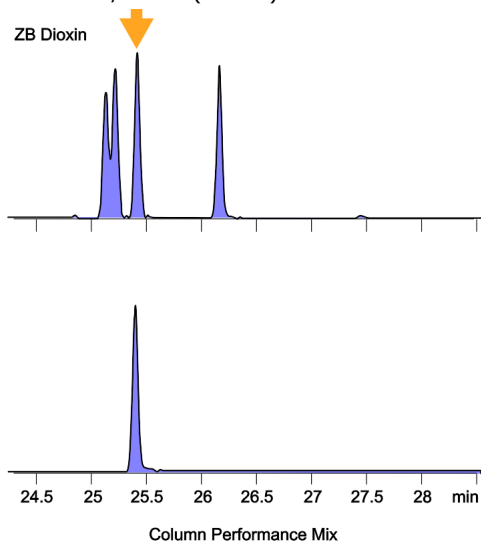


2,3,7,8-TCDD Column Performance and Comprehensive Test Mix on ZB-Dioxin**Column:** Zebron™ ZB-Dioxin, GC Cap. Column 60m, 0.25mm, 0.20um, Ea**Phase:****Dimensions:** 60 meters x 0.25 mm x 0.2 µm**Order No:** 7KG-G045-10**Oven Profile:** 190°C for 3.5 min to 230°C @ 8°C/min 230°C for 15 min to 290°C @ 6.7°C/min 290°C for 5 min to 315°C @ 21°C/min 315°C for 9 min**Carrier Gas:** Constant Flow Helium, 1.3 mL/min**Injection:** Splitless :1 2 µL @ 280°C**Detection:** GC/HRMS (300°C)**Products used in this application:****ANALYTES:**