

Fluconazole Standard solution (mix of 4) USP on Kinetex 5 µm C18 150x4.6mm

Column: Kinetex® 5µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4601-E0

Elution Type: Isocratic

Eluent A: Water:Acetonitrile 80:20 (Pre-mixed)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	15	100


Flow Rate: 0.5 mL/min

Col. Temp.: 40 °C

Detection: LC/UV (DAD, PDA) @ 260.000000000 nm (nanometers) (ambient)



Products used in this application:

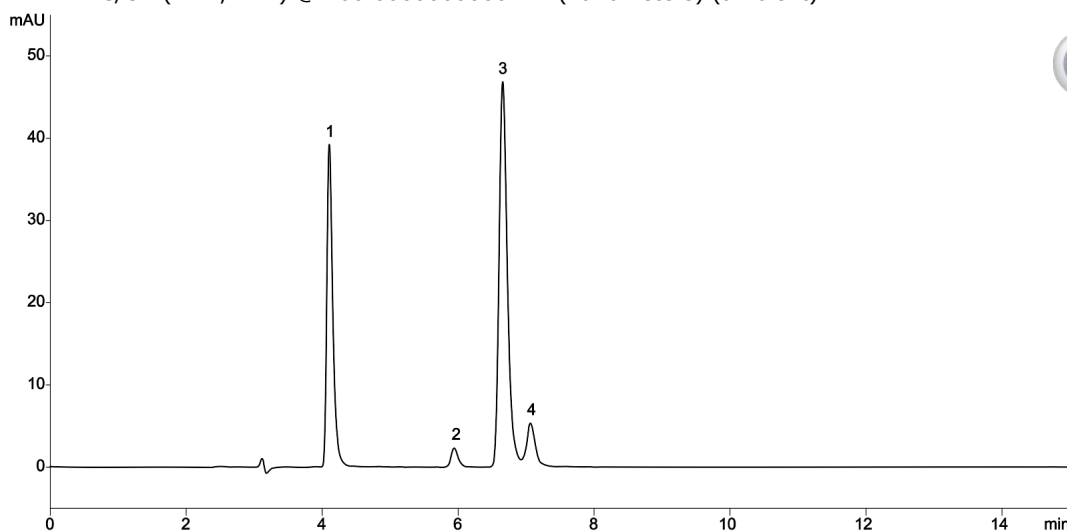

ANALYTES:

- 1 Fluconazole Related Compound A

Retention Time: 4.106 min

- 2 Fluconazole Related Compound B

Retention Time: 5.943 min

- 3 Fluconazole Related Compound C

Retention Time: 6.656 min

- 4 Fluconazole

Retention Time: 7.064 min

