

Fluconazole Standard solution (mix of 4) USP on Kinetex 2.6 µm C18 100x4.6mm

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4462-E0

Elution Type: Isocratic

Eluent A: Water:Acetonitrile 80:20 (Pre-mixed)

Gradient

Step No.	Time (min)	Pct A
1	15	100


Flow Rate: 0.5 mL/min

Col. Temp.: 40 °C

Detection: LC/UV (DAD, PDA) @ 260.000000000 nm (nanometers) (ambient)



Products used in this application:

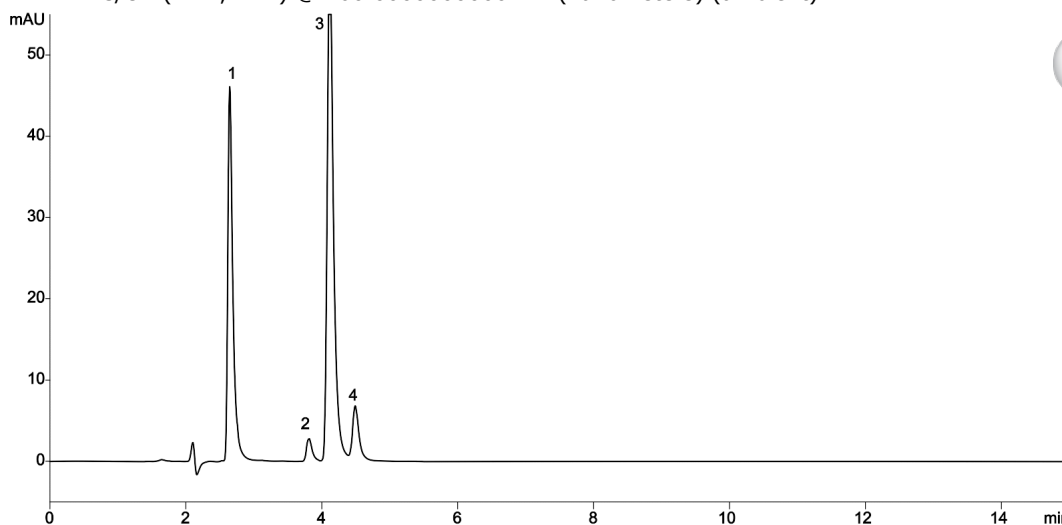

ANALYTES:

1 Fluconazole Related Compound A

Retention Time: 2.649 min

2 Fluconazole Related Compound B

Retention Time: 3.814 min

3 Fluconazole Related Compound C

Retention Time: 4.12 min

4 Fluconazole

Retention Time: 4.493 min

