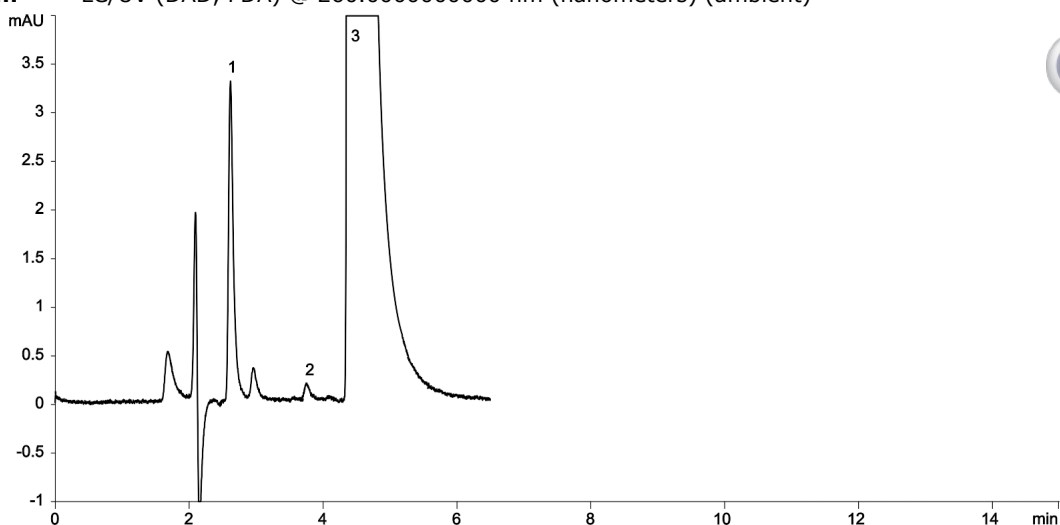


Fluconazole Sample Solution 3 mg/mL on Kinetex 2.6 µm C18 100 x 4.6 mm**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 100 x 4.6 mm, Ea**Dimensions:** 100 x 4.6 mm ID**Order No:** 00D-4462-E0**Elution Type:** Isocratic**Eluent A:** Water:Acetonitrile 80:20 (Pre-mixed)**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 15 100**Flow Rate:** 0.5 mL/min**Col. Temp.:** 40 °C**Detection:** LC/UV (DAD, PDA) @ 260.000000000 nm (nanometers) (ambient)**Products used in this application:****ANALYTES:**

- 1**
- Fluconazole Related Compound A

Retention Time: 2.642 min

- 2**
- Fluconazole Related Compound B

Retention Time: 3.788 min

- 3**
- Fluconazole

Retention Time: 4.432 min

