

# HPLC Application

ID No.: 3265

## Minoxidil - USP Acceptable Modified Method

**Column:** Luna® 5 µm C18(2) 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4252-E0

**Elution Type:** Isocratic

**Eluent A:** Methanol

**Eluent B:** Water

**Eluent C:** Acetic acid with 7mM Docusate sodium, pH 3.0

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C
1	1	0	69.3	29.7	1

**Flow Rate:** 1.5 mL/min

**Col. Temp.:** ambient

**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

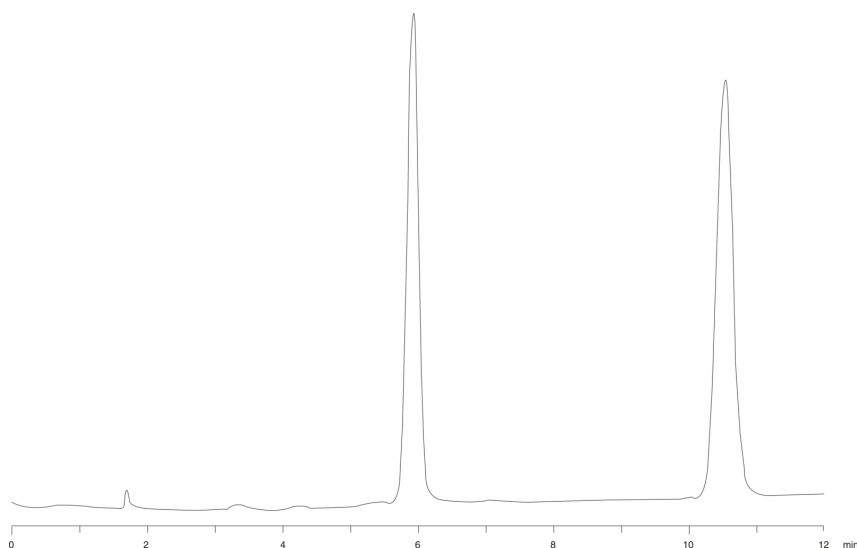

### ANALYTES:

1 Medroxyprogesterone acetate

2 Minoxidil

