

## Glyphosate in Herbicide Roundup

**Column:** Hamilton® 10 µm PRP-X100, LC Column 100 x 4.1 mm, Ea

**Dimensions:** 100 x 4.1 mm ID

**Order No:** CH0-0685

**Elution Type:** Isocratic

**Eluent A:** 40 mM Succinic acid pH 2.9

| Gradient | Step No. | Time (min) | Pct A |
|----------|----------|------------|-------|
| Profile: | 1        | 0          | 100   |

**Flow Rate:** 1.7 mL/min

**Col. Temp.:** ambient

**Detection:** Conductivity @ \_ (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, RP-1 4 x 3.0mm, 10/Pk Part No.: AJ0-5809

- Holder Part No.: KJ0-4282



Products used in this application:

**Hamilton®**

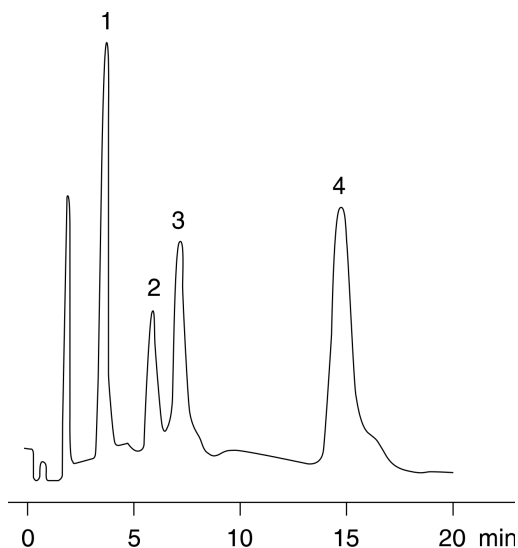

### ANALYTES:

- 1 Fluoride
- 2 Glyphosate
- 3 Phosphate
- 4 Chloride

